

Multi-Application Laser Technology...

With so many patient options in aesthetic medicine today, your practice must be technologically sophisticated to stay competitive. To help meet this need, Focus Medical, a recognized leader in aesthetic lasers, offers the NaturaLase LP -- an incredibly powerful, versatile, reliable, and well designed laser system that creates multiple treatment opportunities from one system. Focus Medical has, once again, raised the bar in Nd:YAG aesthetic laser systems.



Leg & facial vein removal

Vascular Lesions



The NaturaLase LP sets a new standard of care for the treatment of leg veins, facial veins, and vascular lesions.

The NaturaLase LP has variable pulse width, spot size, and pulse rate to effectively treat a wide range of leg veins, facial veins, and vascular lesions. The 1064 wavelength is uniquely suited for leg veins as the penetration depth and low affinity for melanin allow the laser energy to reach much deeper veins than other technologies. Variable pulse widths and spot sizes ensure the most effective treatment settings. Facial vein and vascular lesion treatments are highly effective and reliable when using proper settings.

More energy with a bigger spot size, equals better beam penetration, faster treatments, and better efficacy.

Skin Rejuvenation *Wrinkle Reduction*


Clinical research confirms that the 1064 wavelength is highly effective at collagen remodeling, and thus improving the texture of the skin and the appearance of rhytides.

The 1064 wavelength penetrates deep into the skin to stimulate collagen remodeling *without* the "down time" of ablative skin resurfacing. The NaturaLase LP along side other gentle treatments such as microdermabrasion and glycolic peels, is the wave of the future. Be part of it.

Hair Removal *Power = Efficacy, Speed, and Safety*

The desire to have smooth, hair free bodies drives an enormous billion dollar global industry. Originally attempted with slow, painful electrolysis and thermolysis, lasers have recently taken center stage as the superior technology -- offering safer, faster, and more effective hair removal treatments.

Though shown clinically effective, the downside of laser hair removal can often outweigh the benefits. Not only can lasers produce swelling and unsightly redness, there are many reports of painful, burning treatments that cause blistering and pigmentation problems.

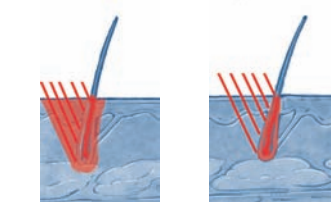

Other Lasers
Shorter wavelength technologies can cause pain, swelling, inflammation or skin damage.

NaturaLase
Energy is focused directly to where it is needed.

Now, the NaturaLase LP offers a breakthrough hair removal solution - a gentle, rapid treatment with minimal patient discomfort that is suitable for *all* skin types.

Long pulse Nd:YAG technology is the most versatile and popular wavelength used in aesthetics. The problem has been speed. With more fluence required to achieve clinical results than other wavelengths, most YAG lasers have been, in comparison, slow. Now, the LP removes the speed barriers to unleash the fastest Nd:YAG hair removal laser ever made.

With 60 watts of power the LP delivers unprecedented speed with fluences and spot sizes that will maximize efficacy and safety. Using a larger spot size increases efficacy and decreases service issues as the number of pulses required per treatment is greatly reduced. Larger spot sizes also require less fluence, which increases safety, especially among darker and tan skinned patients.

The LP is the most versatile, effective, and powerful Nd:YAG available. This translates into satisfied patients and successful practices.

Clinical Results, Unrivaled Safety *Best in Class Nd:YAG Technology*

Highly effective treatment for leg veins. Nd:YAG 1064 is the "Gold Standard" for this application.



Facial Vascular Treatment for a wide variety of facial veins and vascular lesions.



Safe, fast, and gentle, Permanent Hair Reduction for ALL skin types.



Today's patients demand safe, comfortable, and efficacious procedures.

*All From One
Laser...*



Patients get the
results that they
want

Integrated air cooling
system protects skin

Large spot sizes for
consistent and fast
operation

Unmatched Power

Made in USA

*want safer, more
lasers.*

Laser Specifications

Laser Source: Nd:YAG
Delivery System: Fiber Optic Cable
Dimensions: 27"(66 cm) x 17"(43 cm) x 50"(124 cm)
Weight: 110 Lbs. (50 kg)
Power requirements: 110, 200, 220, or 240VAC

NaturaLase LP

Wavelength	1064 nm
Fluence	up to 400 J/cm ²
Spot Size	3 - 15 mm
Pulse Energy	60 J
Pulse Rate	up to 5hz

For More Information Call
(866) MED-LASE
(866) 633-5273 toll free
(203) 730-8885 outside the USA

© 2005 Focus Medical, LLC
NaturaLase® is a trademark of Focus Medical, LLC.
23 Francis J. Clarke Circle
Bethel, CT 06801
www.NaturaLase.com

NaturaLase LP

Multi-Application System

Hair Removal

Leg & Facial
Veins

Wrinkle
Reduction

Vascular
Lesions

NaturaLase LP
Long Pulse Nd:YAG Laser System